THE SMALL MOLECULE SKELETAL SARCOMERE ACTIVATOR, CK-2017357, IS A CALCIUM SENSITIZER THAT BINDS SELECTIVELY TO THE FAST SKELETAL TROPONIN COMPLEX

Raja Kawas, Alan Russell, Alex Muci, Bradley Morgan, Fady Malik, Jim Hartman

Cytokinetics, Inc., South San Francisco, CA, USA

INTRODUCTION

Muscle contraction is driven by the cyclical interactions of myosin thick filaments and actin-containing thin filaments. This process is regulated by both nucleotide-dependent conformational changes in the myosin motor and calcium-dependent conformational changes in the thin filament, controlled by the calcium-sensing protein complex troponin. Muscle types differ in their sarcomere components, particularly the isoforms of myosin and troponin, providing specialization to meet their particular physiological

Direct modulation of muscle contractility at the level of the contractile apparatus is a therapeutic approach with applicability to several diseases. Previous discovery efforts directed at cardiac muscle resulted in the identification of CK-1827452 (now omecamtiv mecarbil), a small molecule direct activator of cardiac myosin that increases cardiac contractility which is currently being studied in Phase II clinical trials in patients with systolic heart failure. Similarly, a small molecule activator of the skeletal sarcomere may have equal utility in increasing muscle function in patient groups where skeletal muscle weakness is a feature.

We have identified a potent troponin activator, CK-2017357, that selectively interacts with the fast skeletal troponin complex, slowing Ca²⁺ dissociation. As a consequence, CK-2017357 sensitizes the sarcomere to calcium, shifting the pCa relationship of both ATPase and force to lower [Ca2+]. We propose a model whereby CK-2017357 sensitizes muscle by stabilizing the Ca2+-bound conformation of troponin, increasing the amount of time thin filaments spend in the active ("open") state.

MATERIALS AND METHODS

Reagents: CK-2017357 was synthesized by the Medicinal Chemistry department at Cytokinetics. Triton X-100 skinned myofibrils were prepared using variations on the published methods of Herrmann et al. (rabbit skeletal) or the working protocol from the lab of Dr. John Solaro (cardiac). Proteins were prepared using variations on the methods of Potter (troponin), Smillie (tropomyosin), Pardee and Spudich (actin), Margossian and Lowey (myosin), and Weeds and Taylor (\$1 fragment digestion).


Myosin isoform analysis was performed using the glycerol SDS-PAGE method of Talmadge and Roy.


ATPase assays were performed at 22°C using a pyruvate kinase/lactate dehydrogenase-coupled enzyme assay system. Skeletal myofibril reactions were performed in a buffer consisting of 12 mM K-PIPES (pH 6.8), 4 mM MgCl₂, 60 mM KCl, 0.5 mM ATP, 0.1 mg/ml bovine serum albumin, 1 mM DTT, 0.5 mM NADH, 1.5 mM phospho-enolpyruvate, 0.6 mM EGTA, and CaCl₂ sufficient to give the desired free calcium concentration. Cardiac myofibril and reconstituted sarcomere assays were performed under similar conditions, except MgCl2 was reduced to 2 mM and KCl was omitted. Reconstituted sarcomere assays contained 0.5 µM S1, 14 μM cardiac actin, 10 μM troponin, and 3 μM tropomyosin.


Isothermal titration calorimetry (ITC) was performed by titrating CK-2017357 (300 µM) into purified rabbit fast skeletal troponin (50 μ M) at 10°C using a VP-ITC calorimeter (Microcal). Both reagents were equilibrated in 12 mM K-PIPES, 100 mM KCl, 250 μM CaCl₂, 5 mM β-mercaptoethanol, and 3% DMSO. Binding parameters were evaluated using a single site binding model.

Calcium release from troponin was measured using the method of Rosenfeld and Taylor (1985), Rabbit fast skeletal troponin (20 µM) or rabbit fast skeletal thin filaments (9 µM troponin, 10 µM tropomyosin, 70 µM actin) were preincubated with either DMSO (2%) or CK-2017357 (40 µM) in DMSO (2%) in PM12 buffer (12 mM K-PIPES pH 6.8, 2 mM MgCl₂, 1 mM DTT) and rapidly mixed with 150 uM Quin-2 in PM12 buffer in a stopped flow apparatus (SF-61DX, TgK Scientific). Quin-2 was excited at 337 nm, with emission measured through a 495 nm long pass filter.

CK-2017357 Selectively Activates Fast Skeletal Myofibrils

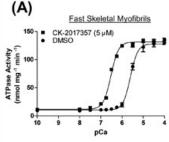

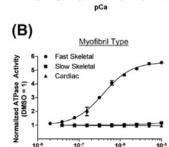

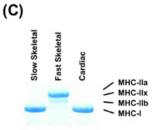

[CK-2017357], mol/L

Figure 1:


CK-2017357 selectively sensitizes the ATPase activity of skinned fast skeletal myofibrils to activation by calcium.

(A) ATPase activity of skinned rabbit fast skeletal myofibrils as a function of free Ca2+ concentration (pCa). The calcium sensitivity shifts from a control pCa50 of 5.7 (DMSO, circles) to 6.8 (5 µM CK-2017357, squares). Error bars indicate the standard deviation of two reactions in a representative experiment.

(B) Sensitization is selective for skeletal myofibrils. CK-2017357 was titrated into skinned myofibrils from fast skeletal (rabbit), slow skeletal (bovine masseter), and cardiac (bovine heart) muscle at approximately the pCa25 for each myofibril type. Error bars indicate the standard deviation of two reactions in a representative experiment.

(C) Myosin isoform analysis showing the myosin composition of the myofibrils used in (B).

CK-2017357 Selectively Activates Thin Filaments Containing Fast Skeletal Troponin

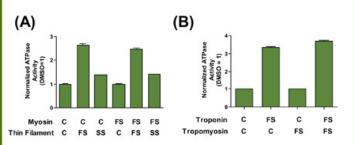


CK-2017357 selectively sensitizes fast skeletal thin filaments to calcium by cting with fast skeletal troponin.

(A) Activation of myosin ATPase is selectively observed with fast skeletal reconstituted thin filaments, regardless of the myosin isoform (fast skeletal or cardiac). The ATPase activity of a "reconstituted sarcomere" system prepared with mixtures of fast skeletal (FS, rabbit), slow skeletal (SS, bovine), and cardiac (C, bovine) components was tested at approximately the pCa25 for each system. All reactions contained cardiac actin. Error bars indicate the standard deviation of two reactions in a representative experiment. CK-2017357 was present at 20 µM.

(B) Activation requires the presence of fast skeletal troponin. The ATPase activity in a reconstituted sarcomere system prepared from combinations of fast skeletal (FS, rabbit) and cardiac (C, bovine) troponin and tropomyosin were tested at approximately the pCa25 for each system. All assays contained cardiac actin and myosin (S1). Error bars indicate the standard deviation of two reactions in a representative experiment.

RESULTS

CK-2017357 Binds to Fast Skeletal Troponin

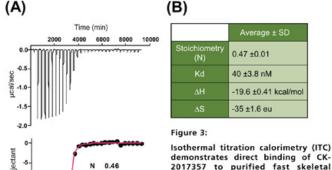

ΔH -19.1 kcal/mol


∆S -33 eu

Molar Ratio

(A) CK-2017357 (300 µM) was titrated into Ca2-saturated rabbit fast skeletal troponin (50 µM) in 3 nmol increments using a VP-ITC calorimeter (Microcal). One repre-sentative experiment is shown, with the results of nonlinear fitting to a single binding site model

(B) The average data from three separate ITC experiments


Potential Binding Model for CK-2017357

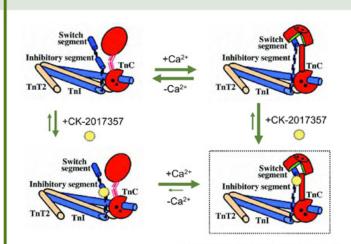


Stabilized Open Conformation

A model for how CK-2017357 could sensitize fast skeletal sarcomeres by stabilizing the open conformation of troponin in the presence of calcium. (Cartoon taken from Vinogradova et al., 2005.)

CK-2017357 Slows Ca2+-Release from Fast Skeletal Troponin and Thin Filaments

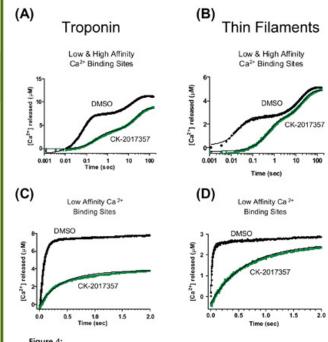

to the 10 µM (troponin) or 4.5 µM (thin filament) final troponin concentration.

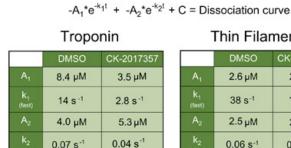



Troponin(Ca2+) + Quin-2 Troponin + Quin-2(Ca2+)

CONCLUSIONS

- reconstituted sarcomere assays containin slow skeletal or cardiac troponin are no activated by CK-2017357.
- 3. CK-2017357 binds directly to purified fast skeletal troponin as demonstrated by isothermal titration calorimetry.
- ponin is slowed by CK-2017357,

This compound has also been shown to rdies of the potential therapeutic uses of ponin activators in diseases where skeletal scle weakness plays a role.

REFERENCES

Davis JP, Norman C, Kobayashi T, Solaro RJ, Swartz DR, and Tikunova SB. Biophysical J. 92:3195–3206(2007).

Herrmann C, Sleep J, Chaussepied P, Travers F, Barman T. A structural and kinetic study on myofibrils prevented from shortening by chemical cross-linking. Biochemistry 32(28):7255-7263(1993).

Potter JD. Preparation of troponin and its subunits, Meth. Enzymo

Rosenfeld SS and Taylor EW. Kinetic studies of calcium binding to

Vinogradova MV, Stone DB, Malanina GG, Karatzaferia C, Cooke R



CK-2017357 slows release of Ca³⁺ from fast skeletal troponin and thin filaments, in particular the regulatory binding sites. Calcium release from fast skeletal troponin (A, C) and thin filaments (B, D) was monitored using the fluorescence enhancement of Quin-2 upon binding to released calcium. Panels A and B show release from both the low (Ca²⁻-specific) and high affinity (Ca²⁻-specific) binding sites. Panels C and D show short time scales where release is primarily from the low-affinity (Ca²⁻-specific) binding sites. Calcium dissociation rates (panel E) were estimated using a biphasic exponential model. Released calcium is expressed as the final concentration after mixing, relative